ITEMNR: 012429
MODELNR: BF6050

HANDLEIDING DIGITALE PERSONENWEEGSCHAAL.....	2
USER MANUAL DIGITAL PERSONAL SCALE	5
BENUTZERANLEITUNG DIGITALE PERSONENWAAGE.....	7
MANUEL D'UTILISATION PÈSE-PERSONNE ÉLECTRONIQUE.....	10

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USER MANUAL DIGITAL PERSONAL SCALE

SPECIFICATIONS

- Large load capacity (180 kg)
- Colour milky white
- Built-in advance step-on technology
- Available units of measurement: kg / lb
- Flat design - only 20mm
- Easy to read LCD display with large digits
- High weighing precision per 100g due to the highly accurate pressure sensor system
- Shows "LO" when the batteries run out
- Shows "OL" when the maximum weight (180 kg) is exceeded
- Switches off automatically after use
- Important: always use the scales on a hard, flat surface, never on carpeting!

Your new electronic personal scales have been designed to accurately indicate your weight gain or loss over a period of time and should give you many years of service under normal use.


GETTING STARTED

1. Please put 2 x AAA Batteries in the battery compartment on the back of the scale.
2. Place the scales on a hard, level surface.
3. Gently step on the scales. Stand on the scales evenly distributing your weight and do not move until your weight reading shown on the display stops changing and gets blocked.

OPERATING INSTRUCTIONS

STEP 1: Always place your scales on a hard, level surface—placing it on carpet or across multiple floor tiles will give you inaccurate results.

STEP 2: Stand on the scales and put your feet evenly on the scale without moving, later after about three seconds your correct weight will be displayed as below picture shows:

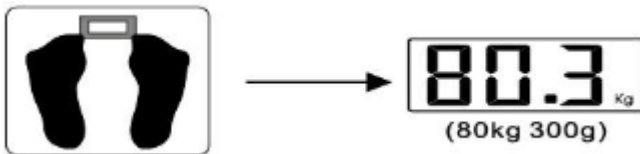
STEP 3: When you step on the scales, it will automatically switch off.

NOTE: When the scale auto-off you just need to operate the scale according to the above mentioned STEP 2 and STEP 3.

WARNING: Keep your scales away from porcelain toilets, sinks, etc. Tempered glass is extremely durable, but can shatter if bumped into porcelain at certain angles.

UNIT BUTTON

The unit button is on the backside of the scale. Press this button to convert unit into KG or LB.

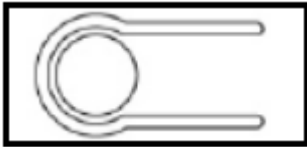BATTERY CHANGE INDICATOR

When the batteries need replacing the following will appear on the display.

Please open the battery compartment in the back of the scales, then take the batteries carefully out of the battery compartment. Then put 2 new AAA Batteries in the battery compartment on the back of the scales.

OVERLOAD INDICATOR

When the scale is overloaded (over 180kg/396LB), the following will appear on the LCD display.

CARE AND ATTENTION

- Place your scales on a hard , flat surface will ensure the greatest accuracy and repeatability.
- Do not stand on the scales with wet feet. The glass on the scales is slippery when wet.
- Always weigh yourself on the same scale on the same floor surface, uneven floors may affect the readings.
- Clean your scales with a damp cloth but do not allow water to get inside the scales.
- Do not use any chemical cleaning agents.
- Do not allow your scales to become saturated with water as this can damage the electronics.
- Remove the batteries if the scales will not be used for over a long period of time.
- Treat your scales with care. It is a precision instrument, do not drop or jump up and down on it.
- These scales are for domestic use only.

STORAGE

Store the product, the manual and accessories in the original packing. This way you always have the necessary information on hand. Store the tool in a dry environment.

DISPOSAL

Do not dispose electric products in the household waste. According to the European standard 2012/19/EU for the disposing of electrical equipment (WEEE), it should be taken to your local waste collection station and offered for recycling.

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INTERBOSCH B.V.